09741

DE

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ab Seite 3

EN

Instruction manual

starting on page 11

FR

Mode d'emploi

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NL

Handleiding

vanaf pagina 27

Z 09741 M DS V4 0318

Inhalt

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
Contents

| | |
|--|----|
| Intended Use _____ | 12 |
| Safety Instructions _____ | 12 |
| Items Supplied and Device Overview _____ | 15 |
| Use _____ | 16 |
| Cleaning and Storage _____ | 17 |
| Troubleshooting _____ | 17 |
| Technical Data _____ | 18 |
| Disposal _____ | 18 |

Explanation of Symbols

Danger symbols: These symbols indicate risks of injury (e. g. from electricity or fire). Read through the associated safety instructions carefully and follow them.

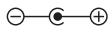Polarity of the coaxial connector (mains adapter)

Protection class II

Protection class III

Circuit symbol for direct current

Energy efficiency class V

May only be used indoors!

Supplementary information

Read operating instructions before use!

Explanation of Signal Words

DANGER

warns of serious injuries and danger to life.

WARNING

warns of **possible** serious injuries and danger to life.

CAUTION

warns of slight to moderate injuries.

NOTE

warns of damage to property.

Dear Customer,

We are delighted that you have chosen this **massage device**.

Before using the device for the first time, please read the instructions through carefully and keep them for future reference. If the device is given to someone else, these instructions should also be handed over with it.

The manufacturer and importer do not accept any liability if the information in these instructions is not complied with.

If you have any questions about the device and about spare parts/accessories, contact the customer service department via our website:

www.service-shopping.de

Intended Use

- This device is intended to be used for loosening the neck and back muscles and delivers relaxation. The device may also be used to massage abdominal, buttock and leg muscles. It is not a replacement for medically required treatment.
- The device is for personal use only and is not intended for commercial applications.
- The device should only be used as described in these instructions. Any other use is deemed to be improper.
- The following are excluded from the warranty: all defects caused by improper handling, damage or unauthorised attempts at repair. This also applies to normal wear and tear.

Safety Instructions

- This device should not be used by people (including children) with impaired sensory or mental faculties or who do not have adequate experience and/or knowledge of such devices, unless they are supervised by a person responsible for their safety, or have received instructions from such a person on how to use the device.
- Children must be supervised to ensure that they do not play with the device.
- The device has a heating function. People who are sensitive to heat must be careful when using the heating function.
- Keep sharp and pointed objects away from the device to ensure that it is not damaged.
- The device must always be disconnected from the mains power when it is left unattended and before cleaning it.

- In the interest of your own safety, you should inspect the device for damage each time before you use it. It may only be used if it has no visible signs of damage!
- If the device, the mains cable or the mains adapter of the device is damaged, it must be replaced by the manufacturer or customer service department or a similarly qualified person (e.g. specialist workshop) in order to avoid any hazards.

Health-Specific Instructions

- **Only use the device if you have no physical illnesses or ailments.**
- The following people may use the device **only after consultation with their doctor**:
 - people who wear a pacemaker or similar medical implants;
 - people with heart disease or who suffer from epilepsy;
 - people who have back complaints or spinal disorders;
 - people who suffer from osteomyelitis, bone inflammation (in particular if the spinal column is damaged) or from osteoporosis;
 - women who are pregnant or have just given birth;
 - people who suffer from circulation problems or thrombosis;
 - people who have undergone an operation for which the treatment is not yet complete;
 - people with physical abnormalities;
 - people who have a malignant tumour.
- If you are unsure whether the device may be used, consult a responsible doctor.
- Do not use the device near to or on injuries, swelling, skin irritations, infected skin, warts, fungal infections, abnormal changes in the skin or body parts with no feeling.
- The device must not be used for longer than 30 minutes continuously. Do not massage one group of muscles for longer than 15 minutes.
- If you experience an unpleasant sensation while you are using the device, stop performing the massage immediately!
- Do not fall asleep while you are using the device!
- Do not use the device before going to bed. The massage has a stimulating effect, which may lead to problems getting to sleep.

DANGER – Danger of Electric Shock from Water!

- Use and store only in closed rooms.
- Do not use the device in rooms with high humidity.
- Never immerse the device, mains cable or mains adapter in water or other liquids and ensure that they cannot fall into water or become wet.
- Should the device fall into water, switch off the power supply immediately. Do not attempt to pull the device out of the water while it is still connected to the mains power!
- Never touch the device, mains cable or mains adapter with wet hands when these components are connected to the mains power.

WARNING – Danger of Electric Shock from Damage and Inappropriate Use!

- Only connect the mains adapter to a plug socket with safety contacts that is properly installed and matches the technical data of the device. The plug socket must also be readily accessible after connection so that the connection to the mains can quickly be isolated.
- Only use suitable extension cables whose technical data is the same as that of the device.
- In the event of damage / faults (e.g. unusual noises, odour or smoke development) and in a thunderstorm, switch the device off immediately and pull the mains adapter out of the plug socket.
- The device must be switched off when inserting or removing the mains adapter. Always pull the mains adapter and never the mains cable to disconnect the device from the mains power.

WARNING – Danger of Injuries!

- **Danger of strangulation!** Make sure that the mains cable is always kept out of the reach of small children and animals.
- The device should only be used and stored out of the reach of children and animals. The mains cable must also be kept out of the reach of children and animals. Otherwise they might use the cable to pull the device down on themselves.
- **Danger of corrosive injury!** (Rechargeable) batteries can be life-threatening if swallowed. Therefore keep (rechargeable) batteries and the remote control out of the reach of young children. If a (rechargeable) battery has been swallowed, you must seek immediate medical attention.
- **Danger of corrosive injury!** If the (rechargeable) battery has leaked, do not allow your skin, eyes and mucous membranes to come into contact with the battery acid. If you touch battery acid, rinse the affected sites with abundant fresh water immediately and contact a doctor as soon as possible.
- Do not insert any fingers or other body parts into the device. Do not open the zip during use.
- **Danger of falling!** Make sure that the connected mains cable does not present a trip hazard.

WARNING – Danger of Fire!

- Do not use the device in rooms containing easily ignitable dust or poisonous and explosive fumes and substances in the air.
- Do not use the device under a blanket or similar item or in the vicinity of combustible material.
- In the event of a fire: Do not extinguish with water! Smother the flames with a fireproof blanket or a suitable fire extinguisher.

NOTICE – Risk of Damage to Material and Property

- Do not lie or stand on the device.
- Protect the device from heat, naked flames, direct sunlight, extreme temperatures, moisture and collisions.
- Do not use any corrosive or abrasive cleaning agents to clean the device. These can damage the surface.
- Keep the device in a dry place at temperatures of between 0 and 40°C.

Items Supplied and Device Overview

PLEASE NOTE!

- **WARNING – Danger of suffocation!** Keep children and animals away from the packaging material.
-
- Unpack all parts and check the items supplied for completeness and transport damage. If you find any damage to the components, do not (!) use them, but contact our customer service department.
 - Remove all possible transport restraints and protective films.

- 1 2 x 2 massage balls, heating and red light element
- 2 Control Panel
- 3 Hand loops

Not shown:

- Mains adapter

Control Panel

The control panel allows you to adjust the following settings:

HEAT (HT) - Red light and heat can be switched off or on again

ON/OFF - The massage device starts the massage and changes direction approx. every minute, the red light and heat function is automatically switched on at the same time

Change of direction (DR) - The massage device massages continuously either clockwise or anticlockwise

Use

PLEASE NOTE!

- No pressure should be exerted on the device as this will slow down or stop the massage and the device could be damaged. Do not stand or lie on the device.
- Do not use the device for longer than 30 minutes without stopping.
- Do not massage one group of muscles for longer than 15 minutes.
- For reasons of hygiene, avoid direct contact of the skin with the device wherever possible. It should not come into contact with skin that has had cream or oil applied to it.
- A tingling sensation/itching during or after the massage is normal due to the circulation-boosting effect.

We recommend that you wait an hour after mealtimes before using the device.

Areas of Application

When this device performs a massage, two massage heads massage the muscles twice. This is intended to relax them and boost circulation. The device is designed in particular for use on the upper back area (in particular the neck). But it can also be used to massage abdominal, buttock and leg muscles.

Neck and shoulders

Back

Abdominal muscles

Buttocks

Thighs

Calves

Operating the Device

1. Connect the mains adapter to the connecting cable on the device.
2. Insert the mains plug into a plug socket that is easily accessible.
3. Place the massage surface of the device onto the part of the body which is to be massaged (see the “Areas of Application” section). For a shoulder or back massage, place the straps of the device around your shoulders and fix them in place by putting your arms through the loops and thus holding the device in the correct position.
4. Start the massage by pressing the **ON/OFF** button.

5. After 15 minutes the device automatically switches off. To switch it off sooner, press the **ON/OFF** button.
6. Pull out the mains plug and store the device in a dry place that is not accessible to children.

Cleaning and Storage

PLEASE NOTE!

- The device may not be connected to the mains and must be completely switched off for cleaning and maintenance!
 - Do not use any caustic or abrasive cleaning agents for cleaning. These can damage the surface.
-
- Wipe down the device just with a damp cloth if necessary. Then wipe it down fully with a dry cloth.
 - Do not wind the mains cable around the device or the mains adapter.
 - Keep the device in its respective packaging in a dry place which is protected from sunlight and is not accessible to children and animals.

Troubleshooting

If the device does not work properly, first check whether you are able to rectify the problem yourself. If the problem cannot be solved with the steps that are listed, contact the customer service department.

Do not attempt to repair a defective electrical device yourself!

| Problem | Possible cause/solution |
|---------------------------|---|
| The device does not work. | <ul style="list-style-type: none"> • Has the ON/OFF button been pressed properly? • Is the mains adapter inserted in the plug socket? • Is the plug socket defective? Try another one. • Check the fuse of the mains connection. • Is the connecting cable connected correctly to the mains cable of the mains adapter? • The device has switched off automatically after operating for 15 minutes. Press the ON / OFF button again. • There is too much weight on the device. |
| No heat is produced. | <ul style="list-style-type: none"> • Is the red light and heat switched on? |
| The device gets too warm. | <ul style="list-style-type: none"> • Was it in operation several times in succession for 15 minutes / too long in each case? Pull the mains adapter out of the plug socket and allow the device to cool down. Always disconnect the device from the mains when it is not in use. |
| The device makes noises. | <ul style="list-style-type: none"> • Operating noises are a normal consequence of the mechanics and not a defect. |

| Problem | Possible cause/solution |
|--------------------------------|--|
| Smoke emerges from the device. | <ul style="list-style-type: none"> Pull the mains adapter out of the plug socket immediately and contact the customer service department. |

Technical Data

| | |
|-------------------|------------|
| Model number: | C37-MP004 |
| Article number: | 09741 |
| Voltage supply: | 12 V |
| Power: | 24 W |
| Timer: | 15 minutes |
| Protection class: | III |

Mains adapter

| | |
|-------------------|--------------------------------|
| Model number: | CP-1225 |
| Input: | 100–240 V ~ 50/60 Hz, max. 1 A |
| Output: | 12 V DC, 2.5 A |
| Protection class: | II |

Disposal

The packaging material is recyclable. Dispose of the packaging in an environmentally friendly manner and take it to a recycling centre.

Dispose of the device in an environmentally friendly manner, it should not be added to general household waste. Dispose of it at a recycling centre for used electrical and electronic devices. More information is available from your local authorities.

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(Calls to German landlines are subject to charges. The cost varies depending on the service provider)

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All rights reserved.
DS Produkte GmbH